

OPTIONS



CB-3003 D ST

PRODUCT SPECIFICATIONS

Operating System Adjustable Hi-Lo, back raise, and foot raise manual by crank

Accessories Side Rail (CSR-001, CSR-002, CSR 003) & Folding Side Rail (FSR-001)

Lock System Double Lock System

Weight (kg) Approx. 90 / **Frame Finishing** Powder Coating

Frame Colour Ivory / **Frame Material** Steel

Dimension (mm) Bed [W 932, L 2141, H 654-994, H/F Board [W 932]

CASTER OPERATION SYSTEM

Independent lock caster (one at front and one at back of the bed) provides stability and secure parking.



FEATURES & PROCEDURES



Head/foot board panel
dismantle



Head/foot board panel
installation

OPERATION SYSTEM

A unique eye catching foldable crank is designed to give high-low bed adjustment.



Crank operation system



Side rail compartment
(available in every
nursing bed)


OPTIONS




CB-3003 D ST

Standard Function

- Operating System : 3 Crank operation for High-Low, Head Rest and Foot Rest positions
- Position : - High – low: Manual operation with Crank
 - Head rest: Manual operation with Crank, tilting from 0 to 80 degrees
 - Foot Rest : Manual operation with Crank, tilting from 0 to 80 degrees
- Primary : - Simple design, no sharp corners (smooth surface)
 - Easy operating mechanism, with Crank System
 - Sturdy and anti-rust
 - Wheels that are easy to move and position, with lock - unlock system
- Lock System : Double Lock System on the Wheels
- Accessories : IV Pole (2 Type Option) and Side Rail / Folding Side rail
 Bed Side Cabinet, Over Bed Table, and Foot Step
 Airmate Cpro Mattress / Urethane

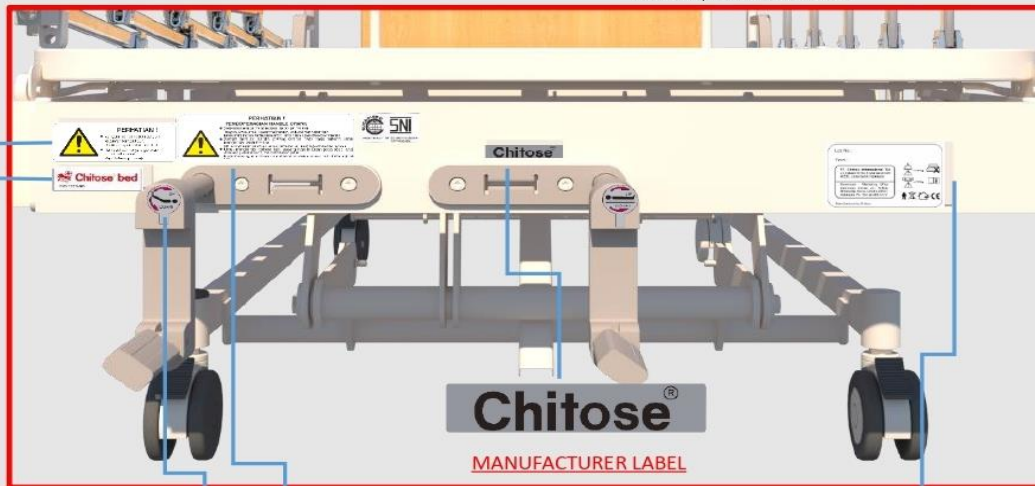

MANUFACTURER LABEL

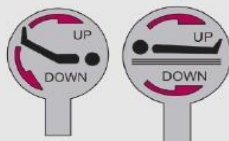

CRANK LABEL



Caution Handle operation



IZIN EDAR KEMENKES



CAUTION LABEL

Lot No :

Type :

PT. Chitose Internasional Tbk
 Jln. Industri III No.5 Utama Cimahi
 40533, Jawa Barat Indonesia

Showroom - Marketing Office
 Distributor Center Jln. H.M.S.
 Mintaredja Baros Cimahi 40521,
 Indonesia Ph. +62 22 206 4777

Manufactured by Chitose

MANUFACTURE LABEL AND LOT NUMBER

- 1 Lot Number and Product Type
- 2 Manufacture Information
- 3 Safe Workload Sign and Maksimum Weight
- 4 Use Appropriate mattress Sign
See Manual Instruction Sign
- 5 Applied Part Sign
Do not Dispose as Household Waste Sign
Indoor Using Only Sign